
## **GI Scope**

1 Leak test before submerging in detergent bath






Check 2.8-5.0 PULL THRU™ Cleaning Brush head against distal end, if too tight or too loose check channel id against scope manufacturer's instructions

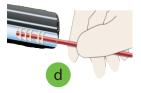




Use 2.8-5.0 PULL THRU™ Cleaning Brush to clean channel x to s

Check for soil on wipers 4 and 5. if there is soil repeat step 3

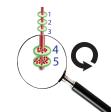


USE 2.8-5.0 PULL THRU™ Cleaning Brush UNTIL THE CHANNEL IS CLEAN



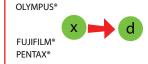

Use 2.8-5.0 PULL THRU™ Cleaning Brush to clean channel x to d

Check for soil on wipers 4 and 5. if there is soil repeat step 4

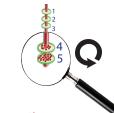

USE 2.8-5.0 PULL THRU™ Cleaning Brush UNTIL THE CHANNEL IS CLEAN





Use 2.8-5.0 PULL THRU™ Cleaning Brush to clean channel y to d

Check for soil on wipers 4 and 5. if there is soil repeat step 5




USE 2.8-5.0 PULL THRU™ Cleaning Brush UNTIL THE CHANNEL IS CLEAN



6



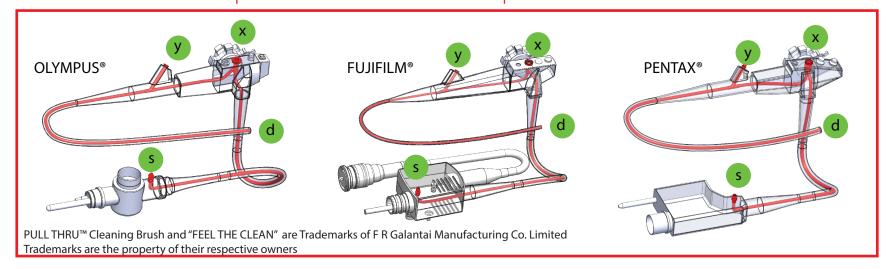

## FORM 55010 Rev 11